SUMMERSHINE SLIM

Lichtdusche
Light Shower
Lámpara terapéutica
Lampe lumineothérapie
Lampada Solare Bianca

10030744

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**Hersteller:**

Zero-Plus International Limited
Room 1004, 10/F., Join-In Hang Sing Centre,
71-75 Container Port Road, Kwai Chung,
New Territories, Hong Kong
E-Mail: sales@zeroplus.hk

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93/42/EWG
2011/65/EU
2017/745/EU

Dear Customer,

Congratulations on purchasing this device. Please read the following instructions carefully and follow them to prevent possible damages. We assume no liability for damage caused by disregard of the instructions and improper use. Scan the QR code to get access to the latest user manual and more product information.

CONTENTS

Technical Data	18
About the Light Shower	19
Signs and Symbols	21
Safety and Performance Standards	21
Safety and Performance Standards (Continued)	23
Safety Instructions	24
Device Overview	26
Commissioning and Operation	26
Cleaning and Storage	27
Troubleshooting	28
Disposal Considerations	28
Declaration of Conformity	29

TECHNICAL DATA

Item number	10030744
Dimensions (W x H x D)	260 x 150 x 17 mm
Weight	450 g
Lighting elements	LEDs
Light intensity	> 10.000 lux in 6 cm
Mains connection	12 V DC 1 A
Operating Environment	0 to +40 °C, and relative humidity from 30 % to 90 %. Atmospheric pressure: 700-1060 hPa
Storage Environment	-10 to +60 °C, and relative humidity of max. 90 %. Atmospheric pressure: 700-1060 hPa
Transport Environment	-10 to +60 °C, and relative humidity of max. 90 %. Atmospheric pressure: 700-1060 hPa
Product Classification	<ol style="list-style-type: none"> 1. Class II equipment 2. No applied part 3. Protection against ingress of water: IP21 4. Not category AP / APG equipment 5. Mode of operation: Continuous operation.
Maximum length of power cord	2000 mm

Included in delivery:

- Light shower
- Power adapter
- Instruction manual

Maximum output of optical radiation for the appliance

Radiance Output	Risk group classified by IEC 60601-2-57:2011	Max. Wert
E_{UVA} : Eye UV-A	Exempt group	$1.0 \times 10^{-4} \text{ W} \cdot \text{m}^{-2}$
ES: Actinic UV skin and eye	Exempt group	$1.0 \times 10^{-4} \text{ W} \cdot \text{m}^{-2}$
EIR: Infrared radiation hazard exposure limits for the eye	Exempt group	$1.0 \times 10^1 \text{ W} \cdot \text{m}^{-2}$
EH: Thermal hazard for the skin	Exempt group	$2.0 \times 10^1 \text{ W} \cdot \text{m}^{-2}$
LB: Blue light	Exempt group	$5.0 \times 10^1 \text{ W} \cdot \text{m}^{-2} \cdot \text{sr}^{-1}$
LR: Retinal thermal	Exempt group	$5.0 \times 10^1 \text{ W} \cdot \text{m}^{-2} \cdot \text{sr}^{-1}$

ABOUT THE LIGHT SHOWER

Notes on operation

Congratulations on the purchase of your light shower. The product is very durable and will provide you with many years of use, provided you handle it correctly. Please keep your receipt, these instructions and the original packaging for reference and storage.

Please read the instructions carefully to get the most benefit from your new device. However, the instructions are not a substitute for advice from a qualified medical advisor.

Intended use

The device is used to treat seasonal depression or winter depression. Daylight increases the energy level of many people and helps them to better adapt to the rhythm of the day and to be more alert. Light showers have been proven to treat the winter blues for several years.

How does the light shower work?

In autumn and winter, the amount of sunlight decreases as the days become shorter. The lack of sunlight causes a problem with certain chemicals in the brain in many people and prevents the hypothalamus from working properly. The lack of light, in return, stimulates the production of the hormone melatonin, which makes you sleepy. The light shower radiates with the strength of sunlight, simulating the lack of light during the dark seasons. At the same time, the dangerous UV rays contained in sunlight are filtered out. The extra light helps the brain reduce and level the production of the hormone melatonin, reducing the symptoms of winter depression.

Symptoms of winter depression

Winter depression can manifest itself through the following symptoms:

- You feel low, listless and fatigued.
- You have trouble concentrating.
- You feel lacking in energy.
- You need more sleep.
- You are less active than usual.

Contraindications

Use is counterproductive and not recommended for the following diseases and symptoms. Always consult your family doctor first before use if you are affected:

- Severe depression.
- Recent eye surgery or diagnoses where bright light should not be allowed to enter the eyes.
- High blood pressure.
- Photosensitive skin or lupus erythematosus.

Possible side effects

The following side effects may occur during treatment with the light shower:

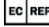
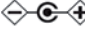
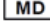
- Headache
- Sleep problems
- Eyestrain
- Nausea

What you can do about it:

- Increase the distance from the light shower.
- Don't use it late at night or just before going to sleep.
- Take a break for a few days until symptoms go away and try again.

SIGNS AND SYMBOLS

The following symbols appear in this manual, on the device or on its accessories:

	Read operating instructions before use		Caution: Consult accompanying documents
	Specifies serial number of the device.		Authorized Representative in the European Community.
	ON/OFF		Class II Equipment
	Date of Manufacture		Manufacturer
	Importeur		DC Polarity
	Medical Device	IP21	Ingress Protection Rating

SAFETY AND PERFORMANCE STANDARDS

This appliance has been designed and manufactured to conform to the highest standards of safety and performance, including electromagnetic compatibility (EMC). This appliance conforms to the applicable requirements of the following:

 This appliance meets the requirement of the Medical Device Directive 93/42/EEC / Medical Device Regulation (EU) 2017/745 for medical devices.

The appliance is complied with the following emission and immunity standard, and test level:

Emisions Test Standard	Compliance
RF emission CISPR 11	Group 1
RF emissions CISPR 11	Class B
Harmonic emissions IEC 61000-3-2	Class A
Voltage fluctuations/flicker emissions IEC 61000-3-3	Comply

Immunity Test standard	IEC 60601-1-2 Test Level	Compliance Level
Electrostatic discharge (ESD) IEC 61000-4-2	± 8 kV contact ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV air	± 8 kV contact ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV air
Electrical fast transient/burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV signal input / output 100 kHz repetition frequency	± 2 kV for power supply lines Not applicable 100 kHz repetition frequency
Surge IEC 61000-4-5	± 0.5 kV, ± 1 kV differential lines ± 0,5 kV, ± 1 kV, ± 1 kV common mode	± 0.5 kV, ± 1 kV differential lines Not Applicable
Voltage dips, interruptions and variations on power supply input lines IEC 61000-4-11	0 % U _T ; 0.5 cycle. At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0 % U _T ; 1 cycle and 70 % U _T ; 25/30 cycles; Single phase: at 0°. 0 % U _T ; 250/300 cycle	0 % U _T ; 0.5 cycle. At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0 % U _T ; 1 cycle and 70 % U _T ; 25/30 cycles; Single phase: at 0°. 0 % U _T ; 250/300 cycle
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m 50 Hz/60 Hz	3 A/m 50 Hz/60 Hz
Conducted RF IEC 61000-4-6	3 V 0.15 MHz – 80 MHz 6 V in ISM and amateur radio bands between 0.15 MHz and 80 Mhz 80 % AM at 1 kHz	3 V 0.15 MHz – 80 MHz 6 V in ISM and amateur radio bands between 0.15 MHz and 80 Mhz 80 % AM at 1 kHz
Radiated RF IEC 61000-4-3	10 V/m 80 MHz – 2.7 GHz 80 % AM at 1 kHz	10 V/m 80 MHz – 2.7 GHz 80 % AM at 1 kHz

Note: U_T is the a.c. mains voltage prior to application of the test level.

SAFETY AND PERFORMANCE STANDARDS (CONTINUED)

Guidance and manufacturer's declaration - electromagnetic immunity							
Radiated RF IEC 61000-4-3 (Test specification for ENCLOSURE PORT IMMUNITY to RF wireless communications equipment	Test Frequency (MHz)	Band (MHz)	Service	Modulation	Maximum power (W)	Distance (m)	IMMUNITY TEST LEVEL (V/m)
	385	380 - 390	TETRA 400	Pulse modulation 18 Hz	1.8	0.3	27
	450	430 - 470	GMRS 460, FRS 460	FM ± 5 kHz deviation 1 kHz sine	2	0.3	28
	710	704 - 787	LTE Band 13, 17	Pulse modulation 217 Hz	0.2	0.3	9
	745						
	780						
	810	800 - 960	GSM 800/900 TETRA 800 iDEN 820 CDMA 850 LTE Band 5	Pulse modulation 18 Hz	2	0.3	28
	870						
	930						
	1720	1700 - 1990	GSM 1800; CDMA 1900; DECT; LTE Band 1, 3, 4, 25; UMTS	Pulse modulation 217 Hz	2	0.3	28
1845							
1970							
2450	2400 - 2570	BT, WLAN, 802.11 b/g/n RFID 2450, LTE Band 7	Pulse modulation 217 Hz	2	0.3	28	
5240	5100 - 5800	WLAN 802.11 a/n	Pulse modulation 217 Hz	0.2	0.3	9	
5240							
5785							

SAFETY INSTRUCTIONS

- Do not use the light shower in damp places.
- Before use, remove all packaging and check the device for damage. If you are unsure, contact customer service.
- The device may only be used by suitable persons.
- Children should be supervised during use. This light shower is not a toy.
- They do not need special training for this device.
- Keep the power cord away from heat and open fire.
- The minimum distance between the light shower and flammable objects must be 1 m.
- The operating voltage must correspond to that indicated on the rating plate.
- After use, disconnect the power plug and allow the light shower to cool down.
- Never touch the light shower with wet hands. Keep them away from moisture.
- Protect the light shower from strong impacts.
- Never pull on the power cord, grasp the plug to pull it.
- If the power cord is damaged, it must be replaced immediately by qualified personnel.
- Never attempt to repair the light shower yourself.
- Disconnecting from the mains can only be done by pulling the mains plug. Always use an easily accessible power outlet.
- Do not use the device near volatile substances, pure oxygen or nitrous gases.
- The device does not require any adjustment.
- The device has no components that can be repaired or maintained.
- The user must not touch the appliance and the person to be treated at the same time.
- The device must not be modified.
- In case of health concerns of any kind, consult your family doctor!

Important notes on electromagnetic compatibility (EMC)

- Use of this appliance adjacent to or stacked with other device should be avoided because it could result in improper operation. If such use is necessary, this appliance and the other device should be observed to verify that they are operating normally. If the appliance does not light on or flickering, it means that the appliance is not operating normally.
- Do not use near active HF surgical equipment and the RF shielded room of an ME system for magnetic resonance imaging, where the intensity of EM disturbances is high.
- Do not change any cable without authorization of the manufacturer. Use of cables other than those specified or provided by the manufacturer of this appliance could result in increased electromagnetic emissions or decreased electromagnetic immunity of this appliance and result in improper operation.
- Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm to any part of this appliance, including power cable. Otherwise, degradation of the performance of this equipment could result.

- If the appliance does not light on or flickering due to the electromagnetic disturbances, do not use it and try to place the appliance in other location.
- For maintaining the basic safety and essential performance with regard to electromagnetic disturbances for the expected service life, precautions shall be taken when the site of use is in the vicinity (e.g. less than 1.5 km away) of AM, FM or TV broadcast antennas.


Important notes on operation

- This device is not intended for commercial or clinical use, but exclusively for personal use in the private household! The use of the daylight shower is NOT a medical treatment.
- If you are taking medication such as painkillers, antihypertensives or antidepressants, you should consult your doctor before using the daylight shower.
- People with retinal diseases should be examined by an ophthalmologist before starting to use the daylight shower, just as diabetics.
- Please do not use in case of eye diseases such as cataract, glaucoma, diseases of the optic nerve in general and inflammation of the vitreous humour.

Maintenance and repairs

- Do not attempt to repair the device yourself. This could result in serious injury. Failure to do so will void the warranty.
- If this device is damaged in any way or does not work, do not continue to use it. Try troubleshooting the problem (see Troubleshooting). If it persists, contact customer service or the dealer, as special tools are required to make repairs. Never attempt to repair the device yourself.
- Do not expose the screen to direct sunlight, spotlights or other sources of heat.
- Do not use alcohol or ammonia-based liquids to clean the screen.

DEVICE OVERVIEW

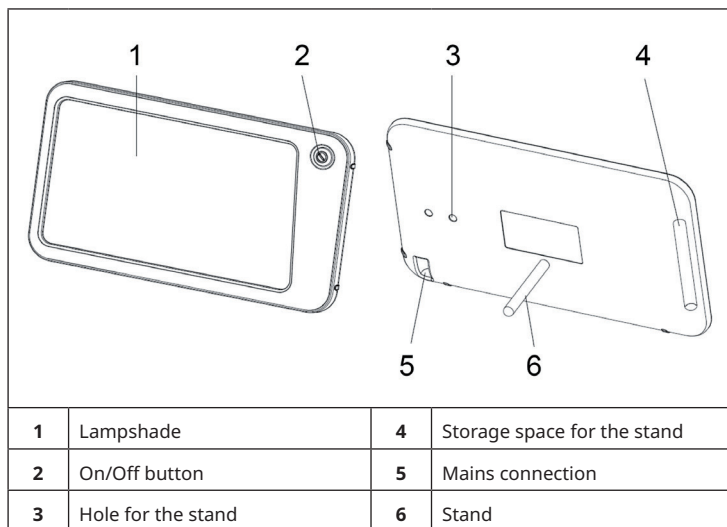COMMISSIONING AND OPERATION

Before operation

1. Remove the stand from the storage compartment.
2. Insert one end of the stand into the hole on the back of the unit.
3. Place the unit on a flat surface about 20 cm from your eyes. The light should hit the side of your face slightly.
4. Plug the device into a 220-240 V outlet. Insert the plug firmly into the outlet and use an outlet that is close to the installation site. Route the cord so that no one will trip over it.

Operation

1. Press a power button to turn on the device.
2. While using, you can read, write, make phone calls or do work. While using the device, try to keep the recommended distance of 20 cm between your eyes and the device. The recommended treatment time is about 2 hours. Important: do not look directly into the light all the time. During the first sessions, you may experience eye pain and headache, which will not occur during the following sessions. In this case, interrupt the application and reduce the duration of use. Do not leave the device unattended during operation.
3. In the darkest season of the year, you should apply the application for at least 7 consecutive days or longer, depending on your personal condition.
4. When you are finished with the application, press the on/off button to turn off the device.

CLEANING AND STORAGE

- Before cleaning, turn the unit off and unplug it from the wall outlet.
- Do not handle the device with wet hands.
- Clean the device with a damp cloth and some washing-up liquid.
- Make sure that no water enters the device.
- Clean the device once a month.
- If you do not use the device for a longer period of time, store it in a dry place. Keep out of the reach of children. Do not leave the device plugged in if you are not going to use it for a long time.

TROUBLESHOOTING

Problem	Possible cause	Solution
No light.	Switched off.	Switch on.
	No power.	Connect to mains.
	No current.	The cable is broken and must be replaced by a specialist company.
	LEDs have reached the end of their service life.	If the device needs to be repaired, contact the manufacturer or a specialist company in your area.
	The LEDs are broken.	

DISPOSAL CONSIDERATIONS

If there is a legal regulation for the disposal of electrical and electronic devices in your country, this symbol on the product or on the packaging indicates that this product must not be disposed of with household waste. Instead, it must be taken to a collection point for the recycling of electrical and electronic equipment. By disposing of it in accordance with the rules, you are protecting the environment and the health of your fellow human beings from negative consequences. For information about the recycling and disposal of this product, please contact your local authority or your household waste disposal service.

DECLARATION OF CONFORMITY

**European Authorised Representative:**

Globalmind Consumer Electronics GmbH
Ernst-Mantius-Str.11,
21029 Hamburg, Germany

**Manufacturer:**

Zero-Plus International Limited
Room 1004, 10/F., Join-In Hang Sing Centre,
71-75 Container Port Road, Kwai Chung,
New Territories, Hong Kong
E-mail: sales@zeroplus.hk

**Imported by:**

Chal-Tec GmbH, Wallstraße 16, 10179 Berlin, Germany.

This product is conform to the following European Directives:

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Chal-Tec GmbH, Wallstraße 16, 10179 Berlino, Germania.